

Contact details

If you are interested in the
KFN-Königsee Femoral nail,
please feel free to contact us



Königsee Implantate GmbH

OT Aschau | Am Sand 4
07426 Allendorf
Germany
export@koenigsee-implantate.de
Fon: + 49 36738 498 - 560
Fax: + 49 36738 498 - 89 560

or your local sales representatives

www.koenigsee-implantate.de



Copyright © Königsee Implantate GmbH/HeBo | Edition 01 | 2020-03

All advantages at a glance

- Nail of 210 mm length; various diameters (10 mm - 13 mm)
- Various CCD angles (125° and 130°)
- Long nail variants
- Colour-coded instrument set for easier handling
- Special supporting screw design with controlled compression path of 15 mm, providing a high degree of primary stability
- Radiolucent light carbon fibre reinforced plastic guide frame KFN with modular technology



Anti-rotation screw, Supporting screw

Sleeve of supporting screw; Locking-, Interlocking-, Compression screw



Indications

Fractures of the femur according to Müller classification AO 31 A1-A3.

Use of the **long KFN interlocking nails** may be considered if the shaft region also has

- Pathological fractures
- Extended subtrochanteric fractures
- Compound fractures (in the proximal femur)
- Ipsilateral trochanteric fractures.

Target group

The KFN is intended for adult patients. It may not be used for children or young people who have not finished growing.



KFN with screws

Features

Product properties - short nails:

- Length: 210 mm
- Cannulated
- Titanium alloy
- anodised surface

Colour coding of nail diameters:

- Ø 10 mm: anodised purple
- Ø 11 mm: anodised green-blue
- Ø 12 mm: anodised blue
- Ø 13 mm: anodised brown

Product variants:

- Item no. 7.1KF.XX/210
· CCD angle: 130°
- Item no. 7.3KF.XX/210
· CCD angle: 125°


Product properties - long nails:

- Ø 10.5 mm
- Cannulated
- Titanium alloy
- anodised goldenyellow

Product variants:

- Item no. 7.1KF10,5L/XXX (left femur)
- Item no. 7.1KF10,5R/XXX (right femur)
- Lengths: 340 mm, 420 mm
- CCD angle: 130°
- Item no. 7.3KF10,5L/XXX (left femur)
- Item no. 7.3KF10,5R/XXX (right femur)
- Lengths: 340 mm, 380 mm
- CCD angle: 125°

Set KFN-Königsee Femoral Nail

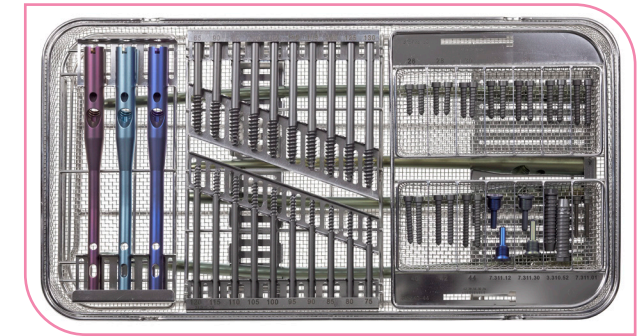


Container OP-Set Implants KFN  Internal hexagon




Item no. 7.1KF.10/210 Item no. 7.1KF.11/210 Item no. 7.1KF.12/210

Container 1 Instruments

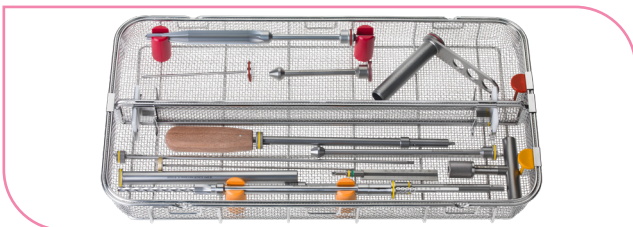

Container 2 Instruments

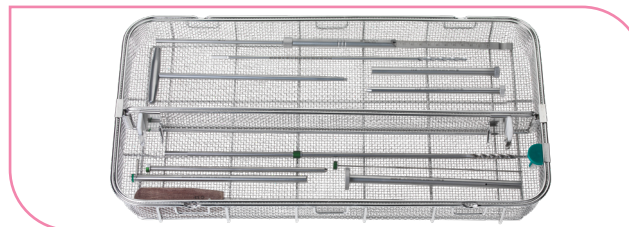

Container 3 Instruments

