

All advantages at a glance

- **Low osseointegration** due to equal distribution of compressive stress over the bones (compared to cables and wires).
- **Studs** on the “bone-facing” side prevent the strip from slipping and also reduce osseointegration.
- **High level of resistance**
- **Metal removal much easier** than with cables or wires – **shortens treatment time**
- **Strip can be retensioned**
- **Universal use** for femur and humerus
- No disruption of **periosteal blood supply**, as proven in extensive studies

Contact details

If you are interested in **LC titanium strip** please feel free to contact us or your local sales representatives.



Königsee Implantate GmbH

OT Aschau | Am Sand 4
07426 Allendorf/Thüringen
Germany

export@koenigsee-implantate.de
+ 49 36738 498 - 560

All implants are available on request as non-sterile packed or sterile packed.



LC TITANIUM STRIP

Low Contact titanium strip with clamp



www.koenigsee-implantate.de

Copyright © Königsee Implantate GmbH
Kr/Mo | Edition 04 | 2023-12



 **Königsee**
Implantate

Target group

The LC titanium strip is intended for adult patients.

LC titanium strip with universal clamping device



Hammer and power cutter




Loop instrument with handle



Pestle



Universal clamping device

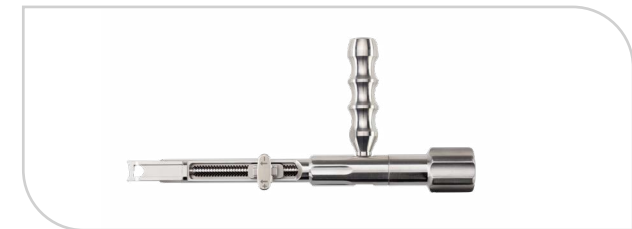

Features

LC titanium strip

Material: pure titanium

The side facing the bone has a studded structure over 70 mm of its length.

There is a recess at one end for hooking the strip into the passer.

The clamp for fixing the strip in place is at the opposite end.


Length 260 mm

Thickness 0.4 mm

Width 8.0 mm

Corresponding instruments and a container to that product is provided in the set LC titanium strip (Set no. 19.519.00).

Set LC titanium strip

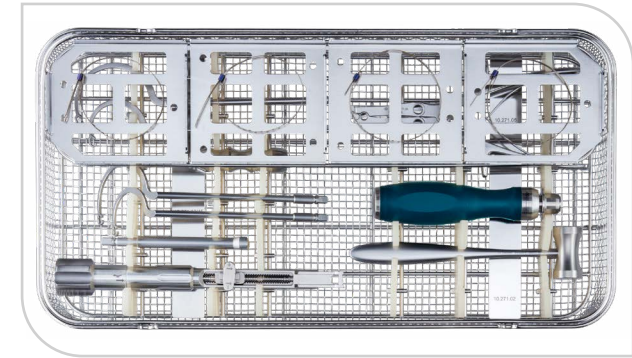

Set LC titanium strip

