

## All advantages at a glance

- Aiming block for the head area for easy installation
- Suitable for deltopectoral approach and deltoid-spitting approach
- Optimal adjustment of plate height using the aiming block
- Temporary fixation of plate and aiming block regardless of screw application
- Improved ability to fix the tubercle to the plate through optimised penetration with surgical needles
- Surgeon can choose between three angle stable screw types for the head holes in the plate head - Ø 3.5 mm cortical / Ø 4.0 mm cortical and Ø 4.0 mm cancellous
- No additional risk to vital ventromedial blood supply
- Access to greater tubercle, frequently subject to dorso-cranial dislocation
- Easy to review the reduction result from cranio-caudal view
- Easy to manipulate fragments using sutures
- Holes are inclined and undercut, making it easy to fix the sutures through the plate
- Can also be used to treat existing (partial) rotator cuff tears

## Contact us

If we have piqued your interest in the **PROXIMALE HUMERUS PLATE EIGHT<sup>2</sup> 4.0/3.5** please contact us or your sales representative.



### Königsee Implantate GmbH

OT Aschau | Am Sand 4  
07426 Allendorf  
Germany  
export@koenigsee-implantate.de  
+ 49 36738 498-560

ALL PRODUCTS ALSO AVAILABLE  
UNDER THE MDR IN NON-STERILE  
AND STERILE.



[www.koenigsee-implantate.de](http://www.koenigsee-implantate.de)


Copyright © Königsee Implantate GmbH  
SiZi | Edition 01 | 2024-08



**Königsee**  
Implantate

## HUMERUS PLATE

Proximal humerus plate Eight<sup>2</sup> 4.0/3.5

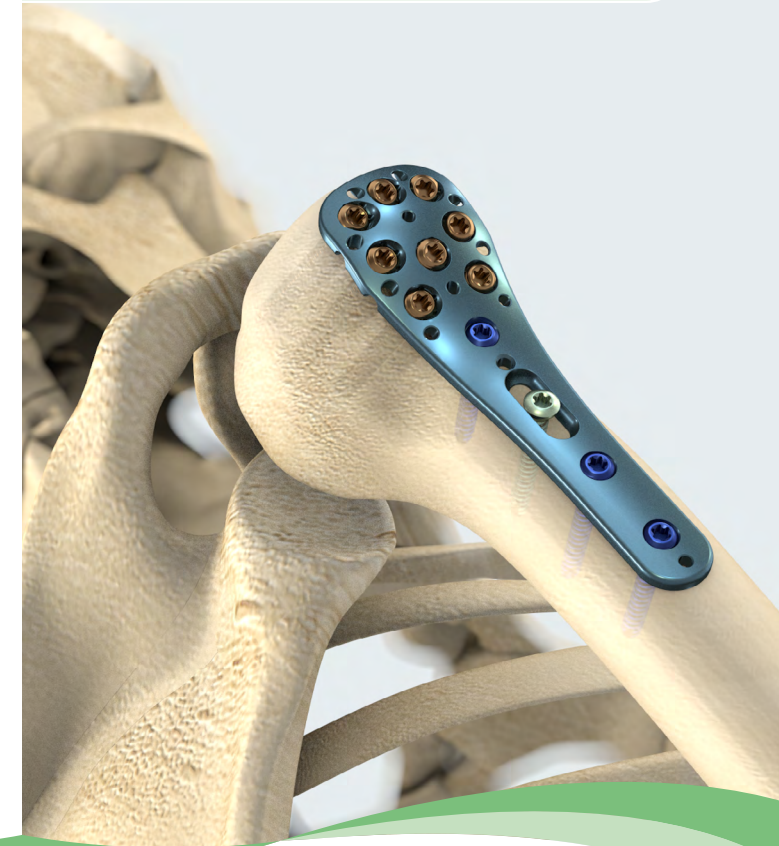

## Screws

### Cortical screw<sup>2</sup> 4.0 as SL15 st



- Ø 2.4 mm
- Self tapping
- Fully threaded
- Available lengths:  
10 - 90 mm
- Anodized surface
- Titanium

### Cortical screw<sup>2</sup> 3.5 as SL15 st



- Ø 3.5 mm
- Self tapping
- Fully threaded
- Available lengths:  
10 - 70 mm
- Anodized surface
- Titanium

SL 15

### Cortical screw 3.5 SL15 st



- Ø 3.5 mm
- Self tapping
- Fully threaded
- Available lengths:  
8 - 140 mm
- Anodized surface
- Titanium

### Cancellous screw<sup>2</sup> 4.0 as SL15 st



- Ø 4.0 mm
- Self tapping
- Fully threaded
- Available lengths:  
10 - 80 mm
- Anodized surface
- Titanium



Can be used optionally  
Aiming block

## Available plate lengths

5.07933.8X  
Proximal humerus plate 4.0/3.5;  
vas; double threaded; titanium

2 hole; length 63 mm

3 hole; length 77 mm

4 hole; length 93 mm

5 hole; length 109 mm

6 hole; length 125 mm


8 hole; length 157 mm

10 hole; length 189 mm


12 hole; length 221 mm

14 hole; length 253 mm

## Sets Proximal Humerus plate EIGHT<sup>2</sup>

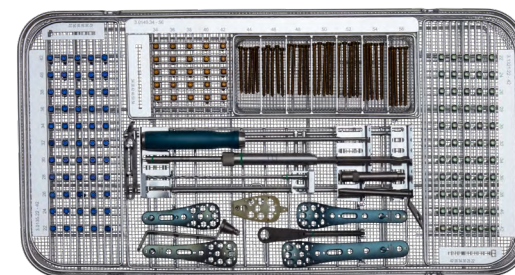

OP-Set Prox.Humerus EIGHT<sup>2</sup> vas Corti

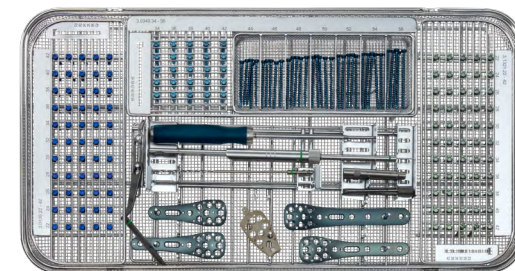

OP-Set Prox.Humerus EIGHT<sup>2</sup> vas



Item no.5.07933.86

